


## DRIVE FULL ELECTRIC LOW HEIGHT BED

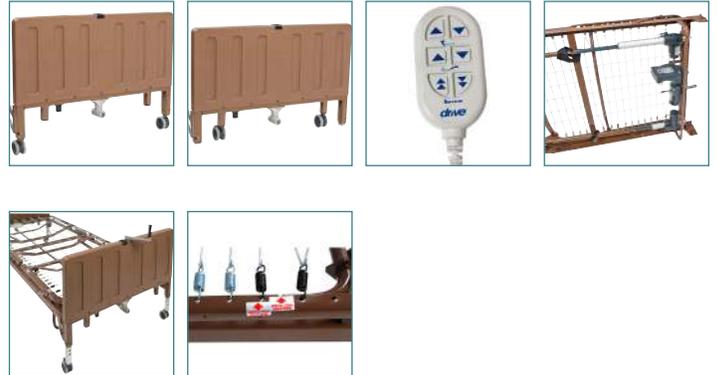

- Low Deck height with flexibility to configure "roll at any height" or "roll at low height only" castor location
- -Roll at any Height: 14.81" (376mm) to 22.83" (580 mm) deck height
- -Roll at Low Height: 11.81 (300 mm) to 19.83" (504mm) deck height
- Quiet, convenient, full electric operation provides smooth position change w/easy-to-use Hand Pendant control of Head/ Foot Section positioning and overall height of the bed frame for maximum comfort and versatility for the patient or care-giver
- Height-adjustment and Head/Foot deck motors come pre-installed on the foot section; no assembly required
- In the event of power failure, top mounted Hi/Lo crank accessibility and manual crank motor override of sleep surface adjustments
- Easy-to-transport and set up, Bed ships in two cartons
- Channel frame construction provides superior strength and reduced weight
- Reinforced frame resists twisting and bending
- Long lasting Zinc-coated spring deck for Patient comfort
- Heavy-duty bed ends with durable HD Polyethylene Molded Bed End covers are attractive, easy-to-maintain
- Headboard is taller than footboard to blend with home décor
- Labels and color-coded springs help ensure accurate installation of side rails for patient safety
- High performance workhorse muscle of the 3-motor power plant is an industry standard
- Half rails and full rails available in packages (PKG) or sold separately
- Motor is UL Approved, Bed is manufactured in ISO9001 & ISO13485 facility
- Satisfies applicable FDA entrapment guidelines
- Extension kit available (model #15005EXTKIT-L) extends length from 80" to 84"
- Patient Weight Capacity 350 lbs. (158kg). Total Weight Capacity 450 lbs. (204kg). Includes occupant, mattress, rails, trapeze, accessories, personal affects, etc.

PRODUCTS	Description
15005LP	Full-Electric Low Height Bed, E0297, 1/cs
15005LP-FR	Full-Electric Low Height Bed with Full Length Side Rails (item #15001ABV), E0266, 1/cs
15005LP-HR	Full-Electric Low Height Bed with Half Rails (item #15201BV), E0266, 1/cs
15005LP-FRPKG	Full-Electric Low Height Bed with Full Length Side Rails (item #15001ABV) and 80" Innerspring Mattress (item #15006), E0265, 1/cs
15005LP-HRPKG	Full-Electric Low Height Bed with Half Rails (item #15201BV) and 80" Innerspring Mattress (Item #15006), E0265, 1/cs
15005LP-FRPKG-TS	Full-Electric Low Height Bed with Full Length Side Rails (item #15001ABV) and 80" Therapeutic Support Mattress (item #15019), E0265, 1/cs
15005LP-HRPKG-TS	Full-Electric Low Height Bed with Half Rails (item #15201BV) and 80" Therapeutic Support Mattress (item #15019), E0265, 1/cs
15005LP-FRPKG-FM	Full-Electric Low Height Bed with Full Length Side Rails (item #15001ABV) and 80" Fiber Mattress (item #15007), E0265, 1/cs

**SPECIFICATIONS**
**Description**

Roll at any Height	14.81" (376mm) to 22.83" (580mm) deck height
Roll at Low Height	11.81" (300mm) to 19.83" (504mm) deck height
Unit Dimensions	36" (.914m) (W) x 88" (2.24m) (L)
Weight	169 lbs.
Carton Shipping Weight	184 lbs.
Weight Limit	Pt Wgt 350 lbs. (158 kg)/Total Wgt 450 lbs. (204 kg) (Includes occupant, mattress, side rails, trapeze, accessories, personal affects, etc.)
Warranty	Original Purchaser. Limited Lifetime on welds, 5 year on frame, motor, pendant, 1 year on all other parts and components
Warranty	During the warranty period, defective items will be repaired or replaced at Drive's option at no charge

**PARTS**
**Description**

15008	Headspring, 1/ea
15013-L	Footspring, 1/ea
15547H-P	Head Board, 1/ea
15015H-FPC	Head Board Cover, 1/ea
15547HPC-1	Head Board Panel, 1/ea
15544HGB	Head Board Gear Assembly, 1/ea
15547F-P	Foot Board, 1/ea
15547FPC-1	Foot Board Panel, 1/ea
15544FGB	Foot Board Gear Assembly (w/o crank), 1/ea
15544C	Hand Crank, 1/ea
15005S	Hi/Low shaft, 1/ea
15005SPN	Adjustment Button, 1/ea
15005HL-L	Hi/Low motor, 1/ea
15005P1029L	Motor Mounting Bracket, 1/ea
15005HLC-L	Hi-Low Motor Cable, 1/ea
15005HLS	Motor Coupling, 1/ea
15015F-FPC	Foot Board Cover, 1/ea
15030P-1001	Caster Receptacle, 1/ea
15044CL	Caster with lock, 1/ea
15044CNL	Caster without lock, 1/ea
15044CSET	Set of 4 Casters for use with Competitor™ Bed (15560), (2 of ea)
15005JBN-L	Junction box, 1/ea
15005JBMB-L	Junction box mounting bracket, 1/ea
15005PC-L	Power cord, 1/ea
15005HCN-L	Hand control, 1/ea
15005MH-L	Head motor, 1/ea
15005P1025-L	Head Motor Tube, 1/ea
15005MF-L	Foot motor, 1/ea
15005PF1025-L	Foot Motor Tube, 1/ea
15005P-1028	Quick Release Spring Knob, 1/ea
15005HFC-L	Head & Foot Motor Cable, 1/ea
15005P1016	Attachment pin & hardware, 1/ea
15005-SPRING	Top Spring, 1/ea
15005-SPRING-1	Side Spring, 1/ea
15005P-1040	Sleep surface hook, 1/bx
15030-FRLINK	Frame Linkage, 1/ea