

Titanium rigid wheelchairs

TiLite<sup>®</sup> ZR | TiLite<sup>®</sup> ZRA



permobil

# TiLite ZR

The 1.25" diameter mono-tube design of the TiLite ZR is the epitome of strength and simplicity. The fixed frame is designed to your specifications and uses aerospace grade titanium for unmatched strength and durability.

Custom designed  
and built to your  
specifications

Fixed or ergo seating  
position is custom  
fabricated to the  
user's specifications

Fixed caster angle saves  
weight and improves the  
chair's responsiveness

Adjustable center  
of gravity allows  
for your perfect  
ride balance



# TiLite ZRA

The classic mono-tube design with a completely adjustable frame. The ZRA is the minimalist look refined.

Custom designed and built to your specifications, ergo seating available

Adjustable center of gravity to find your perfect ride balance

Adjustable seating angle to maximize comfort

1.25" diameter tubing to maximize performance and minimize visual impact

Adjustable caster fork angle to ensure a precision ride



## Tattoo options

- 1 Mosaic
- 2 Flower garden
- 3 Need for speed




## Frame colors

### High gloss paint




### Matte textured paint

### Candy paint

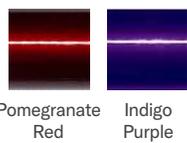

### Metallic paint

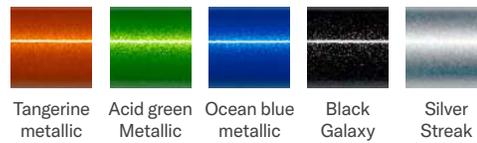

### Matte paint



### Pearlescent paint




### Natural Titanium



## Component colors

### Anodized

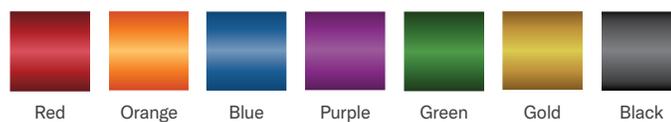

## Specifications

### TiLite ZR

<b>Frame features</b>	Titanium Mono-Tube Rigid, TiFit, Limited adjustability	<b>Center of gravity adjustment</b>	2.5" from 0" – 5" position
<b>Total chair width*</b>	23"	<b>Folding seat back</b>	Adj. height 9.5" – 22" Adj. angle 80° – 101°
<b>Seat width</b>	10" – 20"	<b>Fixed seat back</b>	Fixed height 8" – 20" Fixed angle 80° – 101°
<b>Front seat width (seat taper)</b>	0", 1.5", 2.5"	<b>Transit options</b>	Transport bracket
<b>Seat depth Ergonomic seat</b>	10" – 20" 5" – 9"	<b>Transport weight*</b>	9.3 lb
<b>Front seat height</b>	16" – 21"	<b>Weight limit</b>	265 lb Standard 300 lb Heavy duty
<b>Rear seat height</b>	12.5" – 20"	<b>Warranty</b>	Five years on frame One year on TiLite components
<b>Front frame angle</b>	70°, 75°, 80°, 85°, or 90°		
<b>Camber</b>	0°, 2°, 4°, or 6°		

### TiLite ZRA

<b>Frame features</b>	Titanium Mono-Tube Rigid, TiFit, Full adjustability	<b>Front seat width (seat taper)</b>	0", 1.5", 2.5"
<b>Total chair width*</b>	23"	<b>Folding seat back</b>	Adj. height 9.5" – 22" Adj. angle 80° – 101°
<b>Seat width</b>	10" – 20"	<b>Ergonomic seat</b>	6", 7", 8"
<b>Seat depth</b>	10" – 20"	<b>Transit options</b>	Transit tie-down
<b>Front seat height</b>	16" – 21"	<b>Transport weight*</b>	10.7 lb
<b>Rear seat height</b>	12.5" – 20"	<b>Weight limit</b>	265 lb Standard 300 lb Heavy duty
<b>Front frame angle</b>	70°, 75°, 80°, 85°, or 90°	<b>Warranty</b>	Five years on frame One year on TiLite components
<b>Camber</b>	0°, 2°, 4°, 6°, 8° or 12°		
<b>Center of gravity adjustment</b>	6"		

\*Weights and widths based on 16" x 16" chair with 0° camber, minimum rear wheel spacing and all standard features. Rear wheels are not included in transport weight. For reference, a pair of standard 24" rear wheels weighs 7.2 lb.